


Micro Raman Accessory

Sample-Ref

Reference Sample
for Environment
Change Monitoring

The Sample-Ref accessory permit to apply a **correction** on the data that have been influenced by **environment changes** (room temperature, laser drift, etc.). It is very useful for fine wavenumber shift measurement (stress & strain, graphene, etc.).

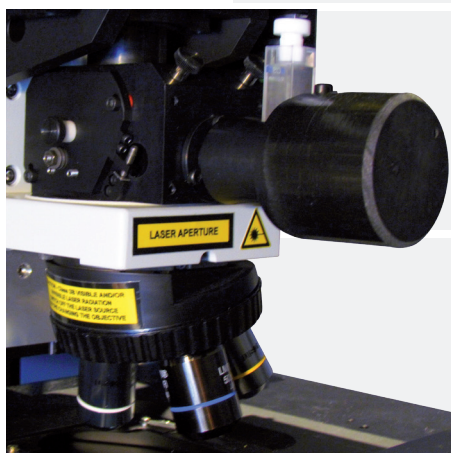

Sample-Ref

The LabRAM HR evolution equipped with this device can exhibit a very high spectral accuracy (typically better than 0.02 cm^{-1}).

- **Spectral drift correction.** A neon lamp spectra can be recorded together with the spectra of the analyzed sample allowing to correct the spectrometer spectral drift. The Raman spectrum of a reference sample can be measured simultaneously to the spectrum of the analyzed sample and permit to correct laser and spectrometer spectral drift.
- **Laser power monitoring.** Measuring a reference sample simultaneously with the sample of interest permits also to monitor and to correct the intensity of variation induced by the laser.
- **Macro Raman.** The same module can also be used as a macro Raman device.


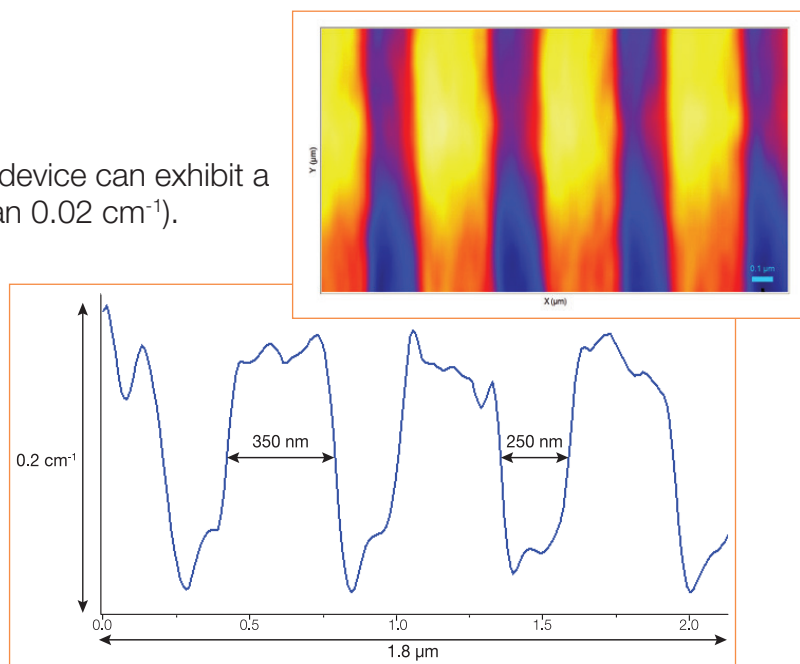


Image and profile of silicon line shift measured on strained silicon

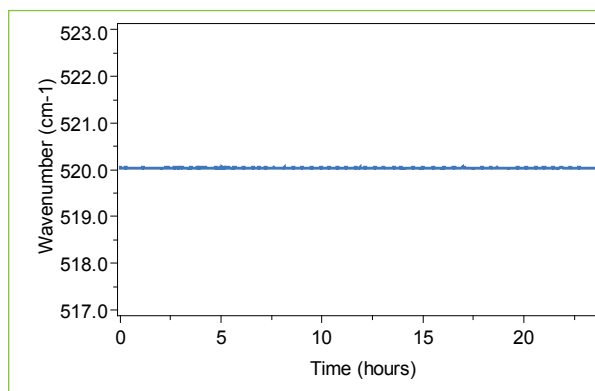

Silicon line position measured over 24 hours using Sample-Ref

This accessory can be easily mounted on the open-space (FSM) upright microscope of the LabRAM HR Evolution.

Find out more at www.horiba.com/raman



HORIBA
Scientific

info.sci@horiba.com

www.horiba.com/scientific

USA: +1 732 494 8660
UK: +44 (0)20 8204 8142
China: +86 (0)21 6289 6060

France: +33 (0)1 69 74 72 00
Italy: +39 2 5760 3050
Brazil: +55 (0)11 5545 1500

Germany: +49 (0)89 4623 17-0
Japan: +81 (0)3 6206 4721
Other: +33 (0)1 69 74 72 00