

Pentra 400 New Generation

The Pentra C400 is an evolution of the successful system, the P400. This evolution ensures the user a platform that is reliable and extremely user friendly, whilst ensuring quality results that is the requirement in such an environment.

Save valuable time

Automatic reagent management

- Intelligent monitoring of reagent supply : the "workbalance" function assesses the stock on board relative to the workload requirements⁽¹⁾
- Use limit date and Calibration validity checks

(1) Workbalance

Real-time workflow follow-up

- Control all your operations from the main screen (system, reagent and sample status)
- Dynamic display of the test results and sampling exceptions with the test review function⁽²⁾
- Multiple alarm sounds available

(2) Test review

Reduced maintenance procedures

- Automatic daily procedures
- Maintenance alerts
- Easy uploading of target values for calibrators and controls

Secure data management

Secured results validation

- Westgard multirule Quality Control procedure with up to 3 controls
- Automatic calibration factor check
- Selective flagging

(3) Quality Control

Reliable traceability in compliance with Accreditation Standards

- Operator login history / Reagents, Calibration and Maintenance logs
- Full patient results traceability : multi-format option for archives review, end of day patient reports
- Full QC program : Up to 12-month control history, Levey-Jennings and Statistics⁽³⁾
- Easy export of archived data (USB key)

Pentra 400

Clinical Chemistry system by HORIBA Medical



Pentra 400

Clinical Chemistry system



HORIBA
Medical



FRANCE +33 (0)4 67 14 15 15 - BENELUX +32 (0)3 281 49 08 - ITALY +39 / 06 51 59 22 1 - SPAIN +34 / 91- 353 30 10 - PORTUGAL +351 / 2 14 72 17 70
UK +44 (0) 1604 542650 - POLAND +48 / 22 6732022 - USA +1 / 949 453 0500 - BRAZIL +55 / 11 5545 1500 - THAILAND +66 / 2 861 59 95
CHINA +86 / 21 3222 1818 - INDIA +91 / 11 4646 5000 - GERMANY AXON LAB AG +49 / 7153 92260 - DISTRIBUTORS NETWORK +33 (0)4 67 14 15 15
HORIBA Medical online : <http://www.horiba.com/medical>



Explore the future

Automotive Test Systems | Process & Environmental | Medical | Semiconductor | Scientific

HORIBA

Explore the future

Automotive Test Systems | Process & Environmental | Medical | Semiconductor | Scientific

HORIBA

Pentra 400

From routine to specific testing

Now, pass your Analysis with flying Colours

- Up to 420 tests per hour with the ISE module
- 55 on-board parameters with back-up availability
- Innovative cassette format
- Integrated work station and validation station
- Touch screen interface
- Real-time follow-up of your workflow



A concentration of technology for high quality results

HIGH-PRECISION OPTICAL SYSTEM

To meet the requirements for highly accurate photometric measurements, HORIBA Medical has chosen the holographic grating technology developed in collaboration with HORIBA Scientific (Jobin Yvon Technology).

- Spectral unit offering 15 different wavelengths (mono or bi-chromatic)
- Optimised reading by continuous sweep
- Highly sensitive measuring (validated holographic technology)
- Space-saving design with maximum efficacy provided by the tungsten lamp

The holographic grating technology reduces light diffraction and provides a high signal / noise ratio which significantly increases the quality of measurement to produce precise and stable test results.



INTEGRATED ISE MODULE

Optional ISE module designed by HORIBA Ltd., Japan.

- 180 tests/hour

- Simultaneous analysis of Sodium, Potassium and Chloride
- Applications on Serum, Plasma and Urine

High quality analytical technology

- Single-use reaction cuvette
To perform the most sensitive tests and all latex assays with high stability results.
- High performance mixer design
For optimal homogenisation in minimum time without cross-contamination.
- Optimum test sequence setting program with Incompatibility function

• High-precision Optical system

• Integrated ISE module

• High quality analytical technology

• An extensive panel of liquid reagents

An extensive panel of liquid reagents

Cassette packaging

- Mono / Bi reagent cassettes
- Ready-to-use
- Extended range
> 52 reagents
> Enzymes / Substrates / Specific Proteins



Easy management of your lab processes

- Random position on board
- Automatic traceability : lot number, expiration date, ...
- Automatic back-up
- Just one calibrator and two controls for global routine analysis



Environment protection

Compact cassette packaging that respects the environment

- No more superfluous packaging
- Saves refrigerator space




Cost control

- Optimised number of tests per cassette
- Limited calibrations
- Re-run analysis according to laboratory defined rules
- Activity measurement thanks to the test counter option⁽¹⁾

Test Counter	Test Counter	Test Counter	Test Counter	Test Counter	Test Counter
1	2	3	4	5	6
7	8	9	10	11	12
13	14	15	16	17	18
19	20	21	22	23	24
25	26	27	28	29	30
31	32	33	34	35	36
37	38	39	40	41	42
43	44	45	46	47	48
49	50	51	52	53	54
55	56	57	58	59	60
61	62	63	64	65	66
67	68	69	70	71	72
73	74	75	76	77	78
79	80	81	82	83	84
85	86	87	88	89	90
91	92	93	94	95	96
97	98	99	100	101	102
103	104	105	106	107	108
109	110	111	112	113	114
115	116	117	118	119	120
121	122	123	124	125	126
127	128	129	130	131	132
133	134	135	136	137	138
139	140	141	142	143	144
145	146	147	148	149	150
151	152	153	154	155	156
157	158	159	160	161	162
163	164	165	166	167	168
169	170	171	172	173	174
175	176	177	178	179	180
181	182	183	184	185	186
187	188	189	190	191	192
193	194	195	196	197	198
199	200	201	202	203	204
205	206	207	208	209	210
211	212	213	214	215	216
217	218	219	220	221	222
223	224	225	226	227	228
229	230	231	232	233	234
235	236	237	238	239	240
241	242	243	244	245	246
247	248	249	250	251	252
253	254	255	256	257	258
259	260	261	262	263	264
265	266	267	268	269	270
271	272	273	274	275	276
277	278	279	280	281	282
283	284	285	286	287	288
289	290	291	292	293	294
295	296	297	298	299	300

(1) Test counter

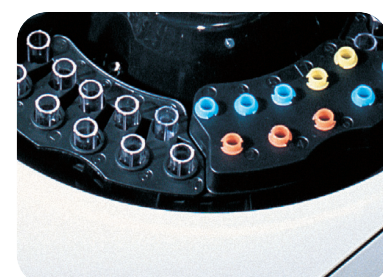


Sample management

- Positive identification
- Multiple sample checks using shock, level and fibrin detection
- Continuous loading
- STAT capability
- Automatic Pre / Post-dilution
- Automatic Re-run



Reagent management

- Barcoded cassettes
- Optimised test packaging
- Reagent test count down with automatic back-up
- Refrigerated zone with up to 44 positions
- Extensive on-board test menu
- User defined channels

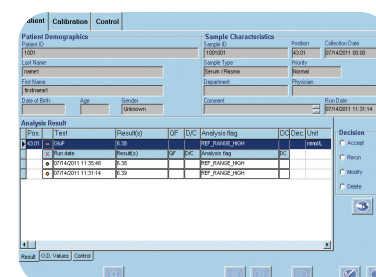

Secure data management

- Secured results validation
- Reliable traceability in compliance with Accreditation Standards



Patient management

- Work list generated by a manual entry or bi-directional connection
- Automatic validation based on user defined rules
- Intelligent delta checking
- Panic values monitoring
- Unlimited patients reports archiving