


## Pentra 400 New Generation

The Pentra C400 is an evolution of the successful system, the P400. This evolution ensures the user a platform that is reliable and extremely user friendly, whilst ensuring quality results that is the requirement in such an environment.

### Save valuable time

#### Automatic reagent management


- Intelligent monitoring of reagent supply : the "workbalance" function assesses the stock on board relative to the workload requirements<sup>(1)</sup>
- Use limit date and Calibration validity checks

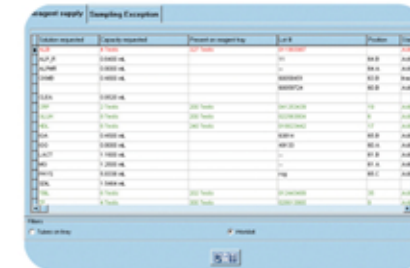

(1) Workbalance

#### Real-time workflow follow-up

- Control all your operations from the main screen (system, reagent and sample status)
- Dynamic display of the test results and sampling exceptions with the test review function<sup>(2)</sup>
- Multiple alarm sounds available

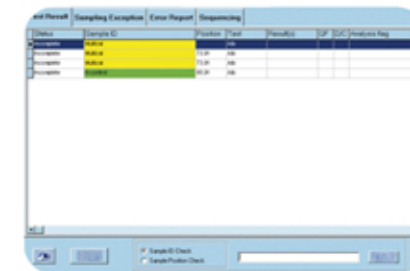

(2) Test review


#### Reduced maintenance procedures

- Automatic daily procedures
- Maintenance alerts
- Easy uploading of target values for calibrators and controls

### Secure data management

#### Secured results validation

- Westgard multirule Quality Control procedure with up to 3 controls
- Automatic calibration factor check
- Selective flagging

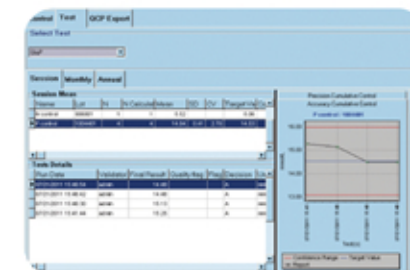

(3) Quality Control

#### Reliable traceability in compliance with Accreditation Standards

- Operator login history / Reagents, Calibration and Maintenance logs
- Full patient results traceability : multi-format option for archives review, end of day patient reports
- Full QC program : Up to 12-month control history, Levey-Jennings and Statistics<sup>(3)</sup>
- Easy export of archived data (USB key)

**Pentra 400**  
Clinical Chemistry system by HORIBA Medical



Subject to technical modifications / HAN751AEN / © 2013 HORIBA Medical - HORIBA ABX SAS - France - RCS Montpellier 328 031 042 - © Shutterstock - James Pangine

**HORIBA**  
Medical

9755 Research Drive, Irvine CA 92618  
949-250-4811 • toll free 1-800-446-7422  
[www.horiba.com/medical](http://www.horiba.com/medical)



**HORIBA**  
Medical

**Pentra 400**  
Clinical Chemistry system





# Pentra 400

## From routine to specific testing

Now, pass your Analysis with flying Colours

- Up to 420 tests per hour with the ISE module
- 55 on-board parameters with back-up availability
- Innovative cassette format
- Integrated work station and validation station
- Touch screen interface
- Real-time follow-up of your workflow

## A concentration of technology for high quality results

### HIGH-PRECISION OPTICAL SYSTEM

To meet the requirements for highly accurate photometric measurements, HORIBA Medical has chosen the holographic grating technology developed in collaboration with HORIBA Scientific (Jobin Yvon Technology).

- Spectral unit offering 15 different wavelengths (mono or bi-chromatic)
- Optimised reading by continuous sweep
- Highly sensitive measuring (validated holographic technology)
- Space-saving design with maximum efficacy provided by the tungsten lamp

The holographic grating technology reduces light diffraction and provides a high signal / noise ratio which significantly increases the quality of measurement to produce precise and stable test results.



### INTEGRATED ISE MODULE

Optional ISE module designed by HORIBA Ltd., Japan.

- 180 tests/hour

- Simultaneous analysis of Sodium, Potassium and Chloride
- Applications on Serum, Plasma and Urine

### High quality analytical technology

- Single-use reaction cuvette  
To perform the most sensitive tests and all latex assays with high stability results.
- High performance mixer design  
For optimal homogenisation in minimum time without cross-contamination.
- Optimum test sequence setting program with Incompatibility function

## An extensive panel of liquid reagents

### Cassette packaging

- Mono / Bi reagent cassettes
- Ready-to-use
- Extended range  
> 52 reagents  
> Enzymes / Substrates / Specific Proteins



### Easy management of your lab processes

- Random position on board
- Automatic traceability : lot number, expiration date, ...
- Automatic back-up
- Just one calibrator and two controls for global routine analysis



### Environment protection

- Compact cassette packaging that respects the environment
- No more superfluous packaging
- Saves refrigerator space




### Cost control

- Optimised number of tests per cassette
- Limited calibrations
- Re-run analysis according to laboratory defined rules
- Activity measurement thanks to the test counter option<sup>(1)</sup>



(1) Test counter

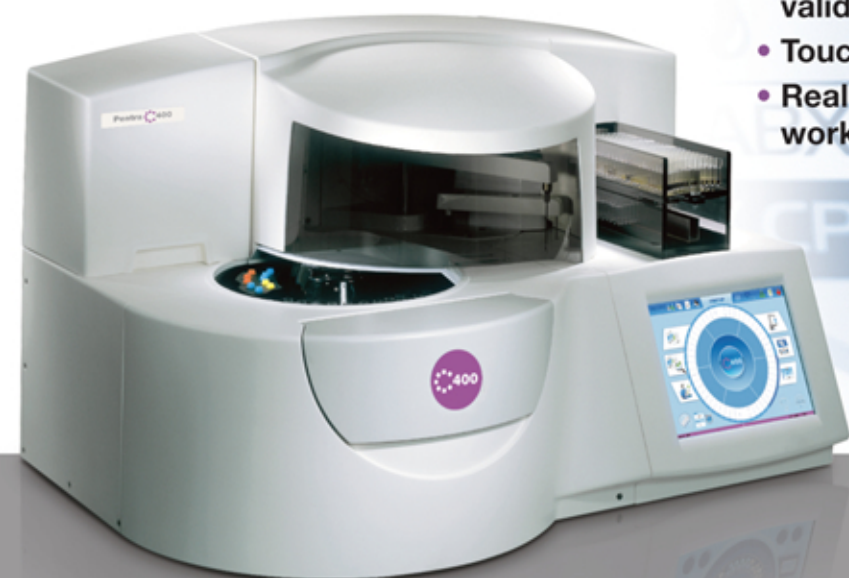

- High-precision Optical system
- Integrated ISE module
- High quality analytical technology
- An extensive panel of liquid reagents



### Sample management

- Positive identification
- Multiple sample checks using shock, level and fibrin detection
- Continuous loading
- STAT capability
- Automatic Pre / Post-dilution
- Automatic Re-run



### Reagent management

- Barcoded cassettes
- Optimised test packaging
- Reagent test count down with automatic back-up
- Refrigerated zone with up to 44 positions
- Extensive on-board test menu
- User defined channels



### Secure data management

- Secured results validation
- Reliable traceability in compliance with Accreditation Standards



### Patient management

- Work list generated by a manual entry or bi-directional connection
- Automatic validation based on user defined rules
- Intelligent delta checking
- Panic values monitoring
- Unlimited patients reports archiving